



STANDARD DEVICE

# DISI-A / Touch

## Automatic disintegration tester

With DISI-A, you control all operation and data analysis safely and easily via computer and PH21 software.


The system fulfils the requirements of FDA 21 CFR Part 11. All operation procedures are password-protected and monitored continuously by the system. All relevant procedures, for example amendments to master data, calibration and much more, are automatically recorded in the log book (audit trail).

The standalone variant DISI-A Touch is easy to operate using the integrated touchscreen. It is possible to upgrade to PC version at any time. The system fulfils all EP/USP requirements.

The device has 4GB of product memory. The standard printout of DISI-A Touch offers an overview of results in table form. You can export data via a USB stick.

	Stations (pieces)	L x W x H (mm)	Weight (kg)
<b>DISI-1A / TOUCH</b>	1	340 x 250 x 750	22
<b>DISI-2A / TOUCH</b>	2	475 x 250 x 750	32
<b>DISI-3A / TOUCH</b>	3	610 x 250 x 750	42
<b>DISI-4A / TOUCH</b>	4	750 x 250 x 750	52

The **DISI-A**, seen here with two measurement stations, is operated via the PH21 computer software

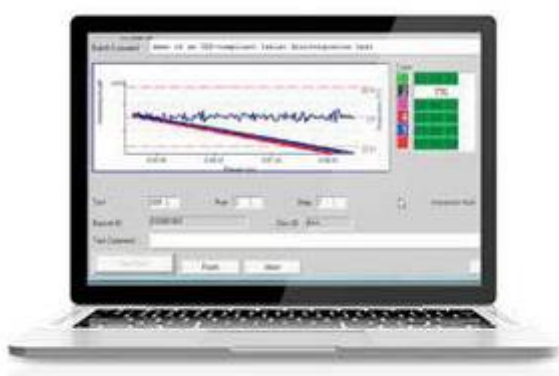

MEASURED PARAMETERS



DISINTEGRATION TIME

ADVANTAGES

- + Touchscreen or computer operation
- + Magnet coupling of the testing baskets
- + Wireless signal transmission
- + Integrated temperature sensors in testing baskets

OPTIONS

- Testing basket B for tablets from Ø 18mm



Special testing basket for large tablets



The **DISI-A Touch** is operated via the touchscreen